



Accessories for EVM series

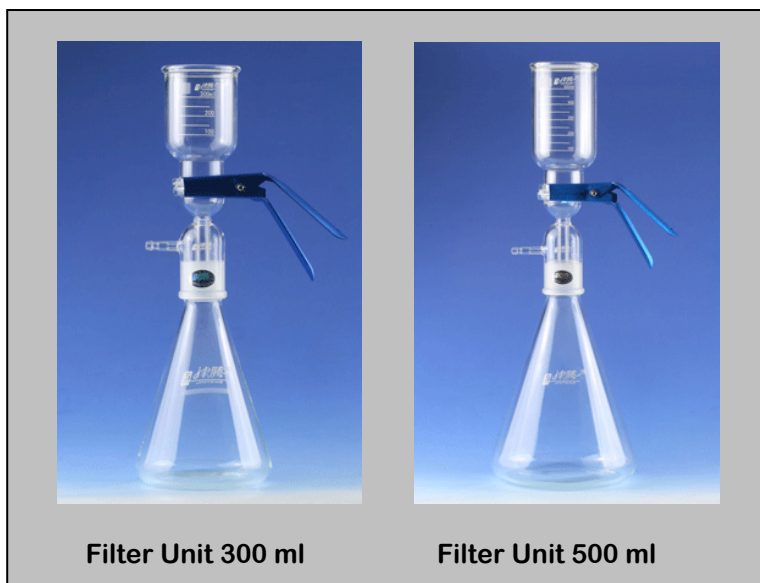

Solvent Filter Unit

The solvent filter is commonly used in the filtration of high-accuracy chromatogram Equipment and has the function of out-gassing to prevent the equipment of blocking caused by particles. But has also many other uses for very fine filtration of liquids. The product is made from high quality extra hard glass with an even wall thickness. It can also be used for sterilisation at high temperature and high pressure. The small unit has a filter cup of 300ml and a receiving bottle of 1000ml – the big unit has a filter cup of 500ml and a receiving bottle of 2000ml. The filter holder has a porous of 10 μ m. The available membrane filters you can find on the next page of this catalog.

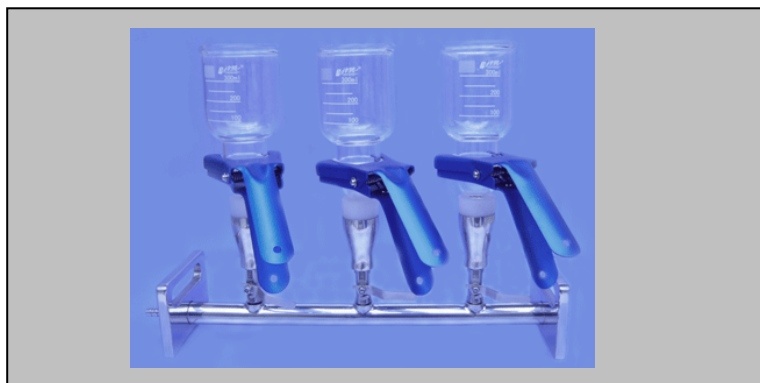

Manifold for 3 solvent Filter Unit

With this manifold it is possible to filter up to 3 filter units on one vacuum pump.

Ordering Information	Solvent Filter Units
Solvent Filter Unit 300ml	629.001
Solvent Filter Unit 500ml	629.002
Manifold Filter Unit	629.003